


SHOULDER ARTHROGRAM – SRS MRI PROTOCOL

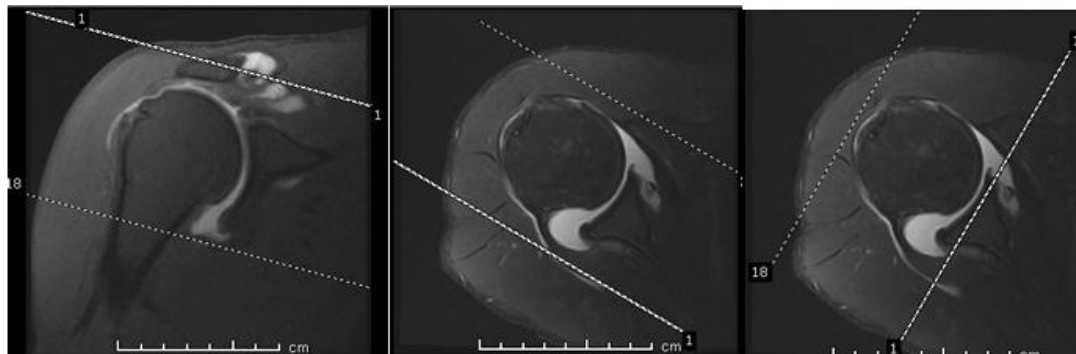

Published: 6/10/2026						Rev: 1			Approval: Dr Tompkins			Magnet/Vendor: Siemens 3.0T XQ			
Seq	Pulse	Dim	Ori	Scan Dir	Phase Enc	FOV (mm)	Slc (mm)	Gap (mm)	TR	TE	Fat Sat	Freq	Phase	Notes	
Localizer	VIBE	2D	Ax	S>>I	A>>P	450	6.0	6.0	7.8	3.69	None	256	138		
SAG & COR AX LOC	VIBE	2D	Sag	R>>L	A>>P	280	6.0	3.0	7.8	3.69	None	256	192		
AXIAL T1 FS	TSE	2D	Ax	S>>I	A>>P	160	3.0	1.0	700	13	FS	320	256		
COR T2 FS BLADE	TSE (BLADE)	2D	Cor	A>>P	A>>P	140	4.0	0.4	3000	43	FS	256	256		
CORONAL T1	TSE	2D	Cor	A>>P	R>>L	160	3.0	0.5	600	11	None	384	346		
CORONAL T1 FS	TSE	2D	Cor	A>>P	R>>L	160	3.0	0.5	700	11	FS	384	346		
SAGITTAL T2 FS	TSE	2D	Sag	R>>L	A>>P	160	3.5	0.5	5050	64	FS	320	256		
SAGITTAL T1 FS	TSE	2D	Sag	R>>L	A>>P	160	3.5	0.5	750	11	FS	320	256		
ABER VIEW LOCALIZER	VIBE	2D	Ax	S>>I	A>>P	450	6.0	6.0	7.8	3.69	None	256	138		
SAG & COR AX LOC (ABER)	VIBE	2D	Sag	R>>L	A>>P	280	6.0	3.0	7.8	3.69	None	256	192		
ABER T1 FS +C	TSE	2D	Sag	R>>L	A>>P	160	3.0	1.0	750	13	FS	320	256		