


WRIST ARTHROGRAM

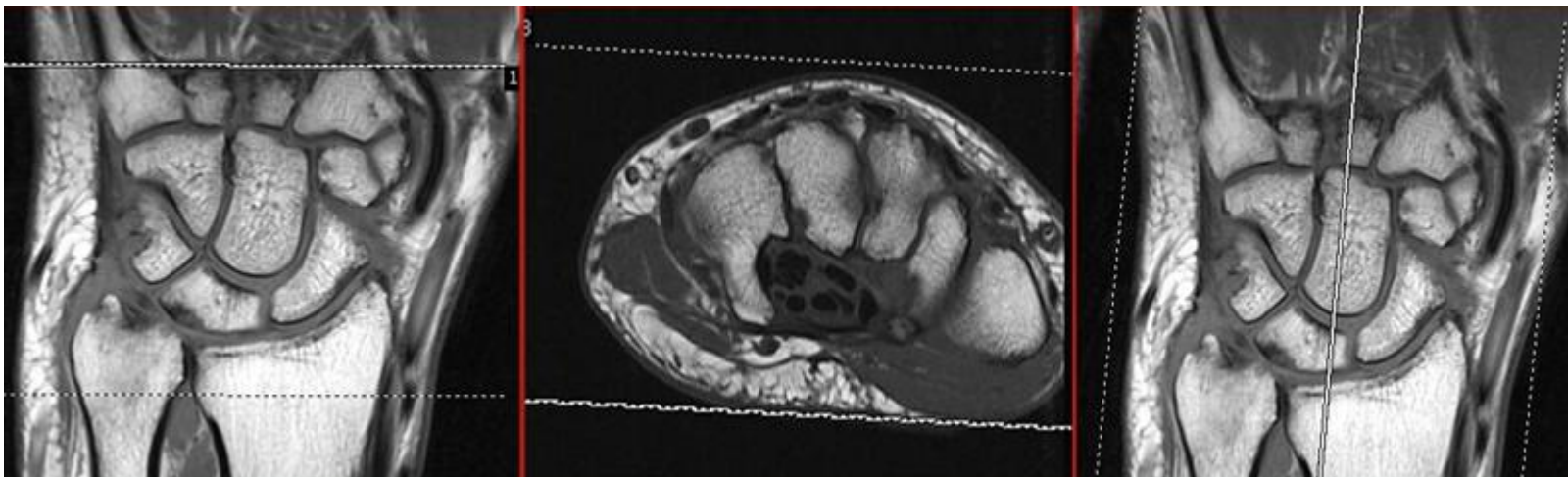

Published: 6/10/2026				Rev: 1			Approval: Dr Tompkins				Magnet/Vendor: Siemens 3.0T XQ			
Seq	Pulse	Dim	Ori	Scan Dir	Phase Enc	FOV (mm)	Slc (mm)	Gap (mm)	TR	TE	Fat Sat	Freq	Phase	Notes
localizer_tra	VIBE	2D	Ax	S>>I	A>>P	400	6.0	3.0	7.8	3.69	None	256	205	
localizer_sag+cor+tra	VIBE	2D	Sag/ Cor/ Ax	R>>L / A>>P / S>>I	A>>P	200	6.0	3.0	8.3	3.69	None	256	205	
AXIAL T1 TSE	TSE	2D	Ax	S>>I	R>>L	90	3.0	0.0	750	13	None	320	240	
AXIAL T1 TSE FS	TSE	2D	Ax	S>>I	R>>L	90	3.0	0.0	700	12	Weak	320	240	
AXIAL T2 FS	TSE	2D	Ax	S>>I	R>>L	90	3.0	0.0	4000	63	Weak	320	208	
CORONAL T1 FS	TSE	2D	Cor	A>>P	R>>L	90	3.0	0.0	600	14	Weak	384	269	
CORONAL T1	TSE	2D	Cor	A>>P	R>>L	90	3.0	0.0	600	14	None	384	269	
CORONAL PD FS	TSE	2D	Cor	A>>P	R>>L	90	3.0	0.0	2500	27	Weak	320	256	
SAGITTAL T2 FS	TSE	2D	Sag	Lat>>Med	H>>F	100	3.0	0.5	4500	64	Weak	320	256	